



Azeer Anesthesia Workstation

HOSPITAL



ANAESTHESIA

Complete range for OR, induction and recovery:



- Anaesthesia devices
- Recovery room
- Accessories

LEON MRI

Specifically designed for use in cardiac cath labs or MRI

Leon mri was specifically developed for use in MRI environments with magnetic field strengths of up to 40 millitesla. It is guaranteed to function with 1.5 Tesla and 3.0 Tesla MRI systems. An LED light that is visible from all sides, including outside of the MR field, clearly displays alarms and notifications in various colours depending on priority.

- The only device in its class that can be equipped with an integrated CO₂ and anaesthesia gas gauge
- Optionally, the device can be connected to a second screen, outside of the MRI range, via a fibre-optic cable
- MRI-compatible touchscreen technology
- Both the performance range and the operation of the Leon mri device match the proven Leon model precisely, without limitations
- The mutual compatibility of the compact circuit systems in the Leon plus, Leon and Leon mri models are another major logistical advantage
- An integrated magnetic field strength monitor verifies the position in relation to the MRI
- The Leon mri model is equipped with a central brake to keep it in position. All four wheels can be locked simultaneously with a single foot lever
- Optional upgrade: Neonatal mode VTG - 3 mlVTG
- Optional upgrade: RS232 data interface via fibre-optic cable

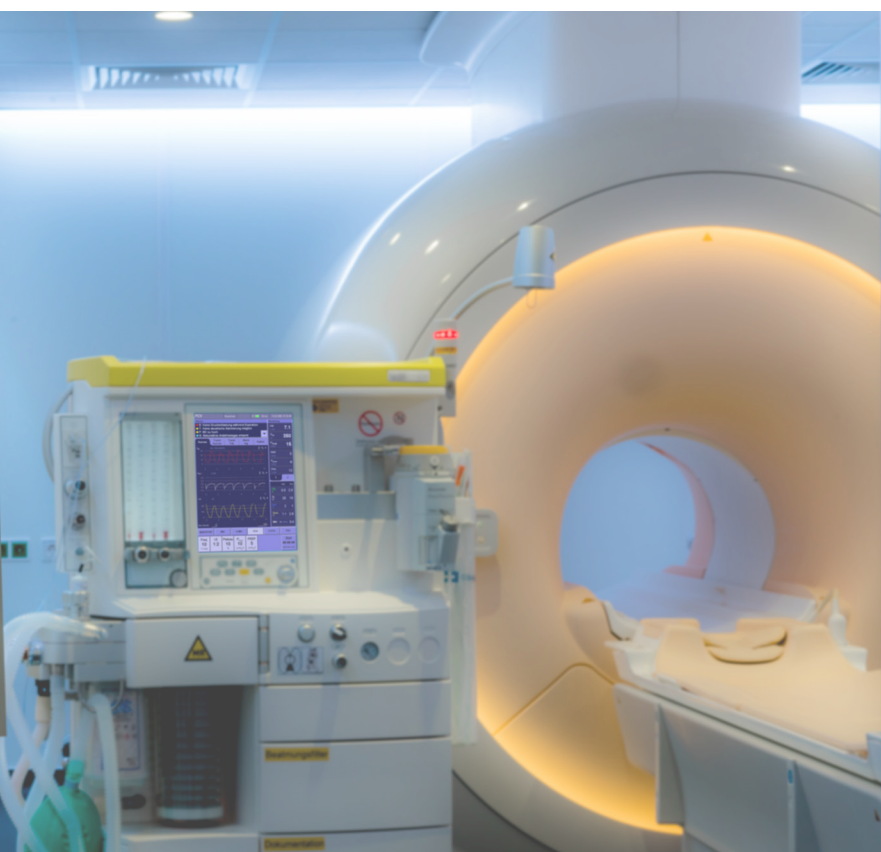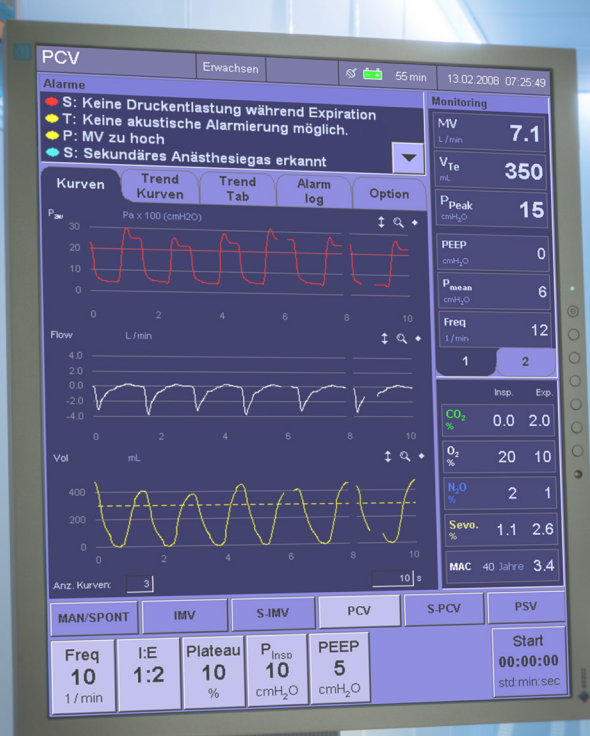

ACCESSORIES

DISPOSABLE ABSORBER

High CO₂ absorption

Leonsorb plus shows highly efficient absorption capacity in low-flow and minimal-flow anaesthesia. Test results have confirmed that the product, with a volume of 1.15 kg, absorbs over 150 litres of CO₂ before a 0.5% CO₂ breakthrough occurs. This property ultimately results in cost savings benefits associated with higher absorption.

The chemical composition of the soda lime used in Leonsorb plus guarantees a high level of CO₂ absorption. A moisture content of approx. 15.5% prevents the soda lime from drying out during operation.



LEONSORB PLUS

- Unique geometry: effective shape for CO₂ absorption
- Reduced response time generates high absorption rates
- Optimum soda lime utilisation
- Minimal dust
- Cost savings due to high absorption capacity and large absorber volume

LEONSORB PREMIUM

- Compatible with low and minimal flow
- Does not generate Compound A
- Unique geometry: effective shape for CO₂ absorption
- Optimum soda lime utilisation
- Minimal dust
- Cost savings due to high absorption capacity and large absorber volume
- Longer lasting colour change than with conventional soda lime





Löwenstein Medical
Arzbacher Straße 80
D-56130 Bad Ems

Tel.: +49 26 03/96 00-0
Fax: +49 26 03/96 00-50
Internet: hul.de



p10319en1802



Images may depict optional accessories that can be ordered separately.

